
# Panasonic x

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# Operating Instructions Instrucciones de funcionamiento

Model No. EW243 Modelo No. EW243

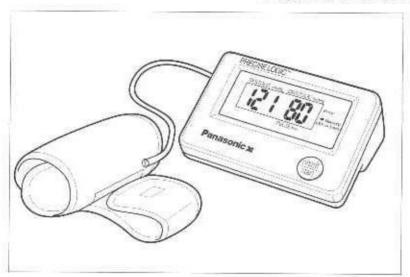Before operating this device, please read these instructions completely.

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Thank you for purchasing the Panasonic Blood Pressure Monitor EW243. The monitor incorporates Precise Logic<sup>™</sup> technology to provide accurate, reliable measurements.

Measuring your own blood pressure is an important way of monitoring your health. High blood pressure (hypertension) is a major health problem which can be treated effectively once detected. Measuring your blood pressure between doctor visits on a regular basis in the comfort of your home, and keeping a record of the measurements using the monitor, will help you monitor any significant changes in your blood pressure. Keeping an accurate record of your blood pressure will help your doctor diagnose and possibly prevent any health problems in the future.

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Your heart acts like a pump, sending blood surging through your blood vessels each time it contracts. Blood pressure is the pressure exerted by blood pumped from the heart on the walls of blood vessels. Systolic pressure is the pressure exerted when the heart contracts and pumps blood into the arteries. Diastolic pressure is the pressure exerted when the heart expands, or relaxes. When you or your doctor take your blood pressure, both your systolic and diastolic pressures are measured. If your blood pressure measurement is 120 over 80 (120/80), for example, your systolic pressure is 120 while your diastolic pressure is 80.

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- Do not confuse self-monitoring with self-diagnosis. Blood pressure
  measurements should only be interpreted by a health professional who is familiar
  with your medical history.
  - No confunda autocomprobación con autodiagnóstico. Las medidas de la presión arterial sólo deberán ser interpretadas por un profesional de la salud que conozca su historial médico.
- If you are taking medication, consult with your physician to determine the most appropriate time to measure your blood pressure. NEVER change a prescribed medication without first consulting with your physician.
  - Si está tomando medicamentos, consulte a su médico para determinar la hora más apropiada para medir su presión arterial. NUNCA cambie un medicamento prescrito sin consultar antes a su médico.
- 3. Blood pressure can vary based on many factors, including age, gender, weight and physical condition. In general, a person's blood pressure is lower during sleep and higher when he or she is active. Blood pressure can change easily in response to physiological changes. The setting in which a person's blood pressure is measured can also affect the results. Having one's blood pressure measured by a healthcare professional in a hospital or clinic can cause nervousness and may result in a temporarily elevated reading. Because blood pressure measurements taken in a clinical setting can vary considerably from those taken at home, a person's blood pressure should be measured not only occasionally in the doctor's office, but also on a regular basis at home. Also, if you find that your blood pressure is lower at home, this is not unusual. To accurately compare with your physician's reading, take your Panasonic blood pressure monitor to your doctor's office and compare readings in this setting.

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- 4. People suffering from cardiac arrhythmia, vascular constriction, liver disorders or diabetes, people with cardiac pacemakers or a weak or irregular pulse, and women who are pregnant should consult their physician before measuring their pressure themselves. Different values may be obtained due to their condition.
  Las personas que sufren de arritmia cardiaca, constricción vascular, enfermedades del hígado o diabetes; las personas con marcapasos cardíacos o pulso débil o irregular; y las
  - Las personas que sutren de arritmia cardiaca, constricción vascular, enfermedades del hígado o diabetes; las personas con marcapasos cardíacos o pulso débil o irregular; y las señoras encintas deberán consultar a sus médicos antes de medir ellas mismas su presión arterial. Puede que se obtengan resultados diferentes según sus condiciones.
- 5. Try to take your blood pressure measurements at the same time and under the same conditions every day.

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#### Blood pressure will be lower than usual:

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 Improper handling of batteries may result in battery explosion or in corrosion from battery leakage. Please observe the followings to ensure proper use of batteries.

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a. Be sure to turn off the power of the device after use.

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d. Be sure to insert batteries with correct polarity, as instructed.

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e. Remove batteries when they are worn out, and dispose of them properly according to all applicable environmental regulations.

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f. Do not disassemble or throw batteries into a fire.

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g. Do not short-circuit batteries.

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h. Do not charge the batteries included with this unit.

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i. If the device requires repressurization to get a proper measurement, there will be a longer measurement time. This may temporarily affect the blood flow in your upper arm, which may feel slightly numb until normal blood flow resumes.

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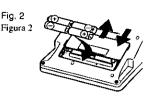
 Lightly press the battery cover on the rear side of the main unit, and push it in the direction of the arrow.

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2. Insert four AA (R6, SUM-3) batteries in the battery compartment, observing the (+) and (-) signs to ensure correct polarity.

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 Approximately 200 measurements can be taken with a new set of Panasonic batteries (at an ambient temperature of 22 °C or 71.6 °F, pressurization of 150 mmHg).

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NOTE:

NOTA:

When the batteries have been almost fully drained, the low battery symbol (Fig.

may not disappear immediately, even after the unit is turned off.

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Either of the above conditions indicates that the batteries are exhausted. You should replace them with four new batteries of the same type.

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NOTE:

NOTA:

Remove the batteries from the unit if it will not be used for an extended period of time (30 days or more).

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### **CAUTION:** ' PRECAUCIÓN:

Do not throw batteries into a fire. Doing so may result in an accident or explosion.

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
# **CAUTIONS:** PRECAUCIONES:

 Do not use the unit for any purpose other than measuring blood pressure. Also, do not use the unit in conjunction with other devices. Improper use may cause a malfunction or personal injury.

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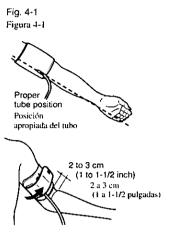
4. Keep the cuff level with your heart (in line with your breast). (Fig. 4-3)

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Fig. 4-3 Figura 4-3
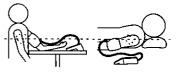
### Measuring Your Blood Pressure Medida de su presión arterial

Make sure the pressure cuff is correctly attached before measuring. The pressure cuff connector at the end of the cuff tube must be firmlg inserted into the pressure cuff connector jack on the main unit.

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#### Relax during measurement.

Relájese durante la medida.

- Press the on/off start button to turn on the unit.
   Presione el botón de encendido/apagado e inicio para encender la unidad.
  - The entire display will to flash for about 2 seconds, and then a zero will be displayed (Fig. 5-1). Inflation will automatically begin immediately thereafter.
     (When the power is turned on and off, you will hear a noise. This is caused by the opening of an air valve and does not indicate a malfunction.)

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
SYSTOLIC moving DIASTOLIC moving formation of the control of the c

2. Measurement begins. (Fig. 5-3)

La medida comienza. (Figura 5-3)

 When inflation has been completed, the measuring indication (♥) will light and the displayed value will gradually decrease.

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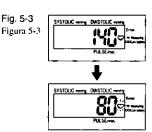NOTE:

NOTA:

The measuring indication ( $\P$ ) will blink in time with the pulse.

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
When measurement has been completed, a completion tone will sound and the air in the cuff will be automatically released.

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4. Blood pressure and pulse values are displayed alternately. (Fig. 5-4)

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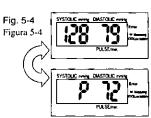(The pulse rate is displayed after converting the number of pulses counted during blood pressure measurement to the pulse rate per minute.)

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#### Error:

#### Error:

If the error indication (♥) next to "Error" (Fig. 5-5) appears on the display during measurement, an error has occurred, e.g., you have moved your hand during measurement. Check to make sure the pressure cuff is positioned and wrapped securely. Check the pressure cuff connection to the main unit, and repeat measurement making sure you remain quiet and still.

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When the correct pulse is not measured, "E"
 (meaning error) appears on pulse rate indication.
 (Fig. 5-6) Noise or an uneven pulse is the probable cause. Please take another measurement.

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5. To take another readings, repeat the procedure from the beginning.

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#### **CAUTION:**

#### **RECAUCIÓN:**

- Do not confuse self-monitoring with self-diagnosis.
   No confunda autocomprobación con autodiagnóstico.
- NEVER change a prescribed medication without first consulting with your physician.

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# Easily Check Your Blood Pressure Readings Against the World Health Organization (WHO) Guidelines

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If your systolic pressure falls in one category but your diastolic pressure in another, your level is classified in the higher of the two categories.

### IMPORTANT: IMPORTANTE:

WHO does not define a low blood pressure area. For reference purposes, systolic
pressure readings of less than 90 mmHg are considered to indicate low blood
pressure regardless of the diastolic pressure.

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Do not be alarmed by temporarily high or low readings because fluctuations in a
person's blood pressure are not uncommon. If possible, measure and record your
blood pressure at the same time every day, and consult your physician if you
have questions or concerns.

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
• If abnormal variations in blood pressure are observed in measurement, please consult your physician.

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The blood flow (pulse wave) which influences blood pressure monitors varies according to the conditions of the blood vessels and the structure of the arm, making measurements by conventional methods difficult in some cases. For more consistent, accurate results, this device employs a Precise Logic™ measurement system that verifies the strength and shape of the pulse wave and then measures the blood pressure.

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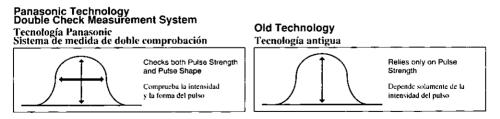By using the Precise Logic<sup>™</sup> measurement system, measurement reliability and accuracy have been improved.

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1. Do not drop the unit. Protect it from sudden jars or shocks.

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2. Do not insert foreign objects into any openings.

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3. Do not attempt to disassemble the unit.

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4. Store this unit out of the reach of children.

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7. If the unit has been stored at temperatures below 0°C (32°F), leave it in a warm place for about 15 minutes before using it. Otherwise the cuff may not inflate properly.

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Do not store the unit in places where there is direct sunlight, high humidity or dust.

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Do not attach the pressure cuff to any part of the body except the arm. NEVER wrap the cuff around the neck.

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 Use a dry, soft cloth to clean the device. If desired, use a cloth lightly dampened with tap water and a mild detergent to clean the device. DO NOT use alcohol, benzene, thinner or other harsh chemicals to clean the device or cuff.

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2. When the unit will not be used for an extended period of time (30 days or more), be sure to remove the batteries. Otherwise, the batteries may leak and damage the unit.

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A replacement pressure cuff for this device can be purchased from an authorized service center.

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If the device still appears to provide unusual or erroneous readings, consult your physician. If the device still does not function properly, contact Panasonic at <u>1-800-338-0552.</u>

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# **Specifications Especificaciones**

 Method of measurement
 : Oscillometric system with Precise Logic™

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Measurement range : 0 - 300 mmHg (Blood Pressure)

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30 - 160 pulsos/minuto (Pulso)

Measurement accuracy : ±3 mmHg (Blood Pressure)

±5% (Pulse rate)

Precisión : ±3 mmHg (Presión arterial)

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Power source : Four AA (R6, UM-3) batteries
Alimentación : Cuatro pilas AA (R6, UM-3)

Operating environment : 10 °C to 40 °C (50 °F to 104 °F) 30% to 85% RH

Ambiente de funcionamiento : 10 °C a 40 °C (50 °F a 104 °F)

30% a 85% de humedad relativa

Storage environment : -20 °C to 60 °C (-4 °F to 140 °F) 10% to 95% RH

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Circumference measurement range : Approximately 20.0 to 40.0 cm (7-3/4" to 15-3/4")

Margen de medida de la circunferencia : 20,0 a 40,0 cm(7-3/4" a 15-3/4") aproximadamente

Dimensions : 69 mm (H)  $\times$ 170 mm (L)  $\times$  115 mm (Dia.)

(2-3/4"(H)×6-3/4"(L)×4-1/2"(Dia.))

Dimensiones :  $69 \text{ mm (altura)} \times 170 \text{ mm (longitud)} \times 115 \text{ mm (diámetro)}$ 

[2-3/4" (altura)  $\times$ 6-3/4" (longitud)  $\times$ 4-1/2" (diámetro)]

Weight : 420g (15 oz) not including batteries

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Accessories : Cuff

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Specifications are subject to change without notice.

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Model EW243 is intended to be used for the oscillometric PRECISE LOGIC™ measurement of systolic and diastolic blood pressure and pulse rate using a pressurized cuff.

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method, within the limits prescribed by the American National Standard, Electronic or automated sphygmomanometers.

# For questions or assistance with your blood pressure monitor, call us at 1-800-338-0552.

Panasonic Personal and Professional
Products Company
A Unit of Matsushita Electric Corporation
of America

One Panasonic Way 4A-3 Secaucus, NJ 07094

© 2000 Matsushita Electric Works, Ltd. Alll Rights Reserved.

No.1 EN, SP (アメリカ) Printed in China